

## USER MANUAL – EN

## IN 24989 Galvanic cosmetic device inSPORTline Hezu




## CONTENTS


SAFETY INSTRUCTIONS	3
PARTS LIST	3
USE	4
RECOMMENDED PROCEDURE	5
SPECIFICATIONS	6
ENVIRONMENT PROTECTION	6
TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS	6

## SAFETY INSTRUCTIONS

- Read the manual before use and keep it for future reference.
- Never repair, modify or disassemble the product.
- Use the device only as instructed in the manual. The device is not intended for any other use.
- Do not store the device in places with high temperature or humidity, such as saunas, bathrooms, etc., there is a risk of electric shock.
- Keep out of reach of children and pets. The device may only be operated by a person over 18 years of age.
- Energy consumption is increased during cooling mode. After using the cooling function, the device must be left to stand for at least 10 minutes.
- Do not cover the ventilation openings.
- Do not use the device if: it emits excessive heat, if the device is damaged, do not use the charging cable if it is damaged.
- Moisturize the skin during use.
- The product is not suitable for: pregnant women, people suffering from epilepsy, allergies to light and heat, people suffering from sensitive skin, dermatological diseases, people with injuries or scars, cancer, or peacemaker etc.

## PARTS LIST

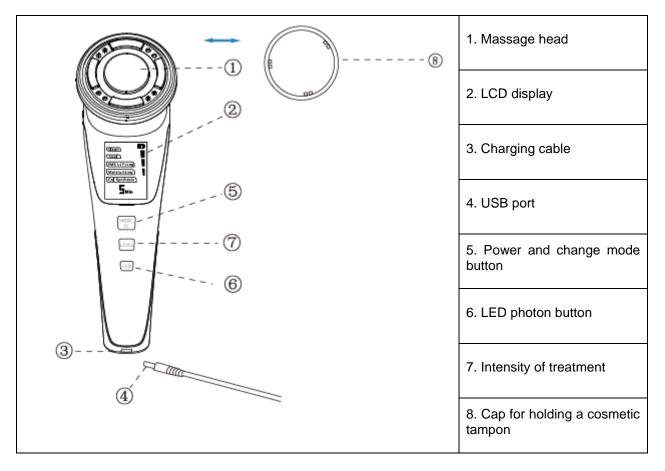

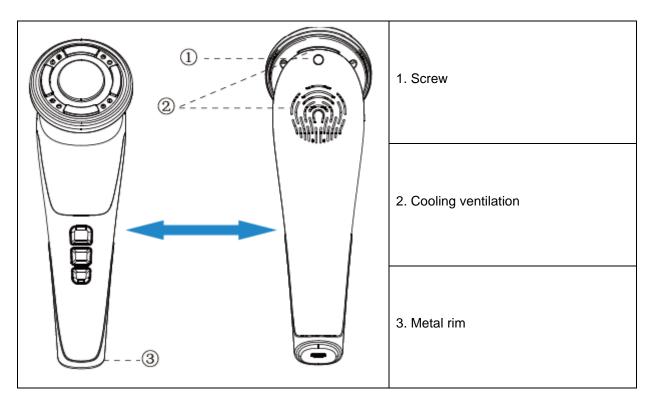

## USE

Make sure the battery is charged before using it for the first time. If necessary, fully charge the battery. Battery charging time is approx. 2.5 hours. Hold down the power button for 5 seconds to start.

- 1. The device starts in vibration and heating mode. You can adjust the temperature with the LEVEL button. You can adjust the LED photonic light using the LED button. For the first function mode, the operating time is set to 5 min. MODE 1.
- 2. Press the power button to select cool mode. The device will vibrate and cooldown. Use the LEVEL button to adjust the temperature. The operating time is 2 minutes, after switching off the fan is in operation for another 2 minutes. MODE 2.
- 3. Use the power button to select EMS lifting mode. The device starts to vibrate and produce an EMS micro current. Use the LEVEL button to adjust the EMS stream production. Operating time is 5 min. MODE 3.
- 4. Use the power button to select moisturizing mode. The device starts vibrating, starts heating and produces an EMS micro current. You can adjust the EMS current strength with the LEVEL button. Operating time is 5 min. MODE 4.
- 5. Press the power button to select the collagen synthesis mode. The device will vibrate, cool and produce an EMS micro current. You can adjust the EMS current strength with the LEVEL button. Operating time is 2 min. MODE 5.

## **RECOMMENDED PROCEDURE**

You can choose between red and blue light for each mode!

## Clean mode:

Cleanse your face from makeup. Start the device, after starting the cleaning function is selected. The cleaning mode opens the pores. The heating function will start. Apply skin cleanser.

#### Cool mode:


Apply a nourishing cream that gets in the open pores. The cooling mode closes the pores and maintains skin moisture. The highest temperature is 42  $^{\circ}$  C, the lowest temperature is 6  $^{\circ}$  C.


#### **EMS lifting mode:**


EMS micro current stimulates muscles and leaves the skin more energetic and elastic. Restores the youth, softness and firmness of the skin.


#### Moisturizing mode:

After cleansing the skin, apply moisturizing cream and use the device to spread the cream over your face. Helps moisturize the skin

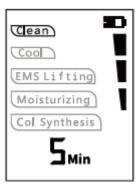

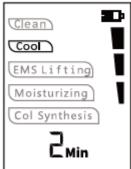


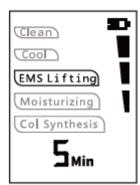

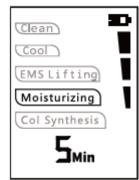

#### Collagen synthesis mode:

After absorbing skin care products, turn on the collagen synthesis mode, the cooling function closes the pores. Increases skin elasticity and cooling closes pores (use skin care products to increase skin elasticity)

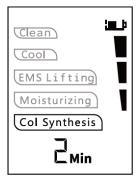

EMS microcurrent stimulates the muscles and skin.

## **SPECIFICATIONS**

Name: Galvanic cosmetic device	Model: IF-1116
Input: DC 5 V	Temperature: 6 - 42 ° C (± 5 ° C)
Vibration: 8000 (± 10) vibrations per minute	Charging time: 2.5 hours
Charging current: 500 mA	Power: max. 8.5 W
Battery capacity: 1000 mAh	Standby time: 45 days
Weight: 145 g	Work time: 5 min or 2 min
Dimensions: 175 x 55 x 67 cm	Package dimensions: 215 x 150 x 80 mm

## **ENVIRONMENT PROTECTION**

After the product lifespan expired or if the possible repairing is uneconomic, dispose it according to the local laws and environmentally friendly in the nearest scrapyard.

By proper disposal you will protect the environment and natural sources. Moreover, you can help protect human health. If you are not sure in correct disposing, ask local authorities to avoid law violation or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

## **TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS**

## **General Conditions of Warranty and Definition of Terms**

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Strakonická street 1151/2c, Prague 150 00, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

"The Buyer who is the End Customer" or simply the "End Customer" is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

"The Buyer who is not the End Customer" is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless

otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

#### Warranty Conditions

#### Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

#### Batteries

6-month battery warranty – we guarantee that battery's nominal capacity does not fall below 70% of its total capacity within 6 months of the product's sale.

#### The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

#### Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such

agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

# CZ SEVEN SPORT s.r.o.

Registered Office: Headquaters: Warranty & Service:	Strakonická 1151/2c, Praha 5, 150 00, ČR Dělnická 957, Vítkov, 749 01 Čermenská 486, Vítkov 749 01
CRN: VAT ID:	26847264 CZ26847264
Phone: E-mail:	+420 556 300 970 eshop@insportline.cz reklamace@insportline.cz servis@insportline.cz
Web:	www.inSPORTline.cz
SK inSPORTline s.r.o.	

Headquaters, warranty & service center: Električná 6471, Trenčín 911 01, SK

CRN:	36311723
VAT ID:	SK2020177082
Phone: E-mail:	+421(0)326 526 701 objednavky@insportline.sk reklamacie@insportline.sk servis@insportline.sk

Web: www.inSPORTline.sk