



USER MANUAL – EN

**IN 26386 Compression massage system inSPORTline
Compisimo**

**IN 31291 Compression Massage System inSPORTline
Compisimo with Unisex Compression Pants - 6
Overlapping Chambers**



CONTENTS

- SAFETY INSTRUCTIONS..... 3
- PRODUCT DESCRIPTION 3
- LED CONTROL PANEL 4
- MESSAGE TYPE 5
- USE 5
- TROUBLESHOOTING 6
- ENVIRONMENT PROTECTION 6
- TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS..... 7

SevenSport s.r.o. reserves the right to make any changes and improvements to its product without prior notice. Visit our website www.insportline.eu where you will find the latest version of the manual.

Before use, it is necessary to unscrew the safety screws, which is intended for safe transportation.

SAFETY INSTRUCTIONS


- Please read the manual before use and keep it for future reference.
- Do not use in a room where there are flammable, combustible or explosive gases.
- Do not use if you suffer from: acute inflammation of the skin, arrhythmia, erythema disease, venous thrombosis, pulmonary edema, acute venous thrombosis, unstable hypertension, people with a pacemaker, varicose veins, have metal implants or other electronic devices.
- Possible side effects: uncomfortable if the pressure is high, redness, itching, pain if the pressure is too high.
- Not suitable for people with limited motor or mental abilities.
- Do not repair or modify the device yourself.
- If you feel discomfort or pain, stop using the product immediately.
- Do not use continuously for more than 20 minutes. Consult your doctor about use.
- Keep out of reach of children and pets.
- Keep away from water and moisture.
- The power cord should be plugged into an easily accessible outlet.
- Power supply 220-240V.
- Keep your hands dry when handling the power cord.
- Protect the power cord from damage.
- Keep away from heat sources or other metals.
- Do not use solvents or other aggressive cleaning agents for cleaning.
- Never use the product on bare skin.
- This is not a medical tool.
- Remote control batteries: CR2025

PRODUCT DESCRIPTION

	
<p>Main unit</p>	<p>Compression for legs</p>
	


Compression for hands – optional	Waist compression – optional
	inSPORTline Unisex Compression Recovery Pants - 6 Chambers

LED CONTROL PANEL

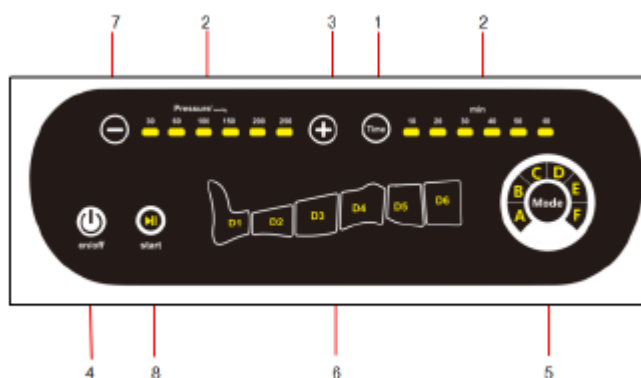

1. Timer 10, 20, 30, 40, 50, 60 min	5. Program type A - F
2. Intensity indicator, 8 levels (0-250 mmHg)	6. Massage position
3. Increase intensity	7. Decrease intensity
4. Power supply	8. Start / pause

The time indicator shows the remaining time of the massage 

By pressing the  button, you can adjust the training time from 10 to 60 min. Each press adjusts the value by 10 min.

Indicator and intensity setting

Use the  and  buttons to adjust the intensity by 10 mmHg, the original set value is 120 mmHg.

Massage position

The default setting is massage using all 6 chambers. Press the button to start/stop the given chamber.

Power switch

Press the button to switch to standby mode, press again to turn off.


Start / pause

In standby mode, the button is green.

If the massage is running, the button is blue.

Press the button to pause or start the massage.

Massage type

Press the button to select the type of massage A - F .

MASSAGE TYPE

Type A: relaxing massage.


Type B: helps relieve the pain of varicose veins.

Type C: deep massage for muscle regeneration.

Type D: helps lose weight in the upper and lower limbs and the waist.

Type E: helps blood circulation.

Type F: combination of massage types A+C+D.

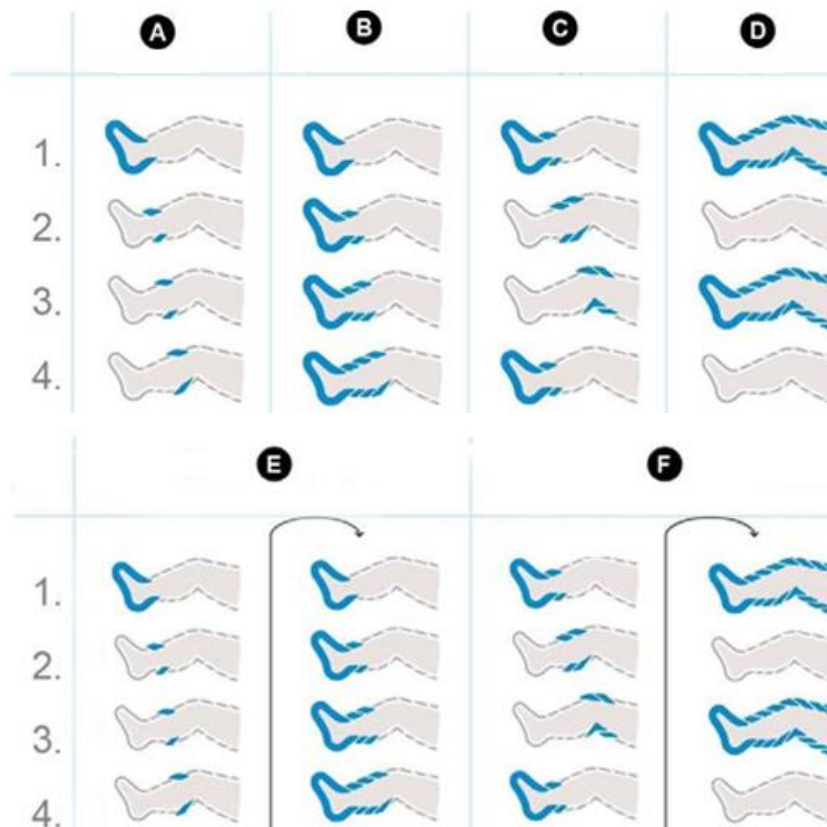

USE

The transport screw must be unscrewed before use. Open the cover and unscrew the screw.



- Place the device on a flat and firm surface.
- Plug in the power and connect the tubing to the device and the given compression sleeve.
- Compression for legs:
- To avoid uncomfortable feeling, make sure the compression sleeve is folded.
- Wear a thin layer of clothing.

- Close the full-length zipper. Connect the hose to the sleeve and to the device. The device should always face the display towards you.
- Turn on the device, massage time, intensity, number of chambers and type of massage. Press the start button to start the massage. (standby mode - green color, in operation - blue color)

Step 1 Plug into a 220-240V power supply.	Step 2 Check the hose connection, if connected skip Step 3.
Step 3 Connect the hose to the compression unit and the device.	Step 4 Make sure all hoses are connected correctly.
Step 5 Nasadte si kompresní rukáv. Varování: Hadici si nedávejte kolem krku.	Step 6 Press the start button.
Step 7 Adjust the intensity using the + and - buttons. Note: In the beginning, we recommend using a lower pressure of 120-150 mmHg.	Step 8 Set the timer. Note: Range 10 – 60 min. Original value 10 min. If the timer reaches 0, the device will turn off.
Step 9 Select a mode.	Step 10 Select the number of chambers. Note: Select the number of massage chambers by pressing the button.
Step 11 If the machine is in operation, press the start/stop button to stop the massage. Note: Do not remove the compression sleeve while the device is paused. The device must be switched off before removal.	Step 12 When finished, disconnect the hose.

TROUBLESHOOTING

The device will not start	<p>Check the power cord connection.</p> <p>Check if the device is running.</p> <p>Check the current voltage.</p> <p>The device automatically turned off after the timer went off</p>
Unnatural sounds (hissing, etc.)	<p>Check the connection and condition of the hose.</p> <p>Check the condition of the connectors</p>

ENVIRONMENT PROTECTION

After the product lifespan expired or if the possible repairing is uneconomic, dispose it according to the local laws and environmentally friendly in the nearest scrapyard.

By proper disposal you will protect the environment and natural sources. Moreover, you can help protect human health. If you are not sure in correct disposing, ask local authorities to avoid law violation or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS

General Conditions of Warranty and Definition of Terms

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Kutnohorská 531, Štěrboholy Retail Park, 109 00, Praha 10, ČR, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

“The Buyer who is the End Customer” or simply the “End Customer” is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

“The Buyer who is not the End Customer” is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

Warranty Conditions

Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

Batteries

6-month battery warranty – we guarantee that battery's nominal capacity does not fall below 70% of its total capacity within 6 months of the product's sale.

The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

CZ

SEVEN SPORT s.r.o.

Registered Office: Kutnohorská 531, Štěřboholy Retail Park,
109 00, Praha 10, ČR

Headquarters: Dělnická 957, Vítkov, 749 01
Warranty & Service: Čermenská 486, Vítkov 749 01

CRN: 26847264
VAT ID: CZ26847264

Phone: +420 556 300 970
E-mail: eshop@insportline.cz
reklamace@insportline.cz
servis@insportline.cz

Web: www.inSPORTline.cz

SK

Stores inSPORTline SK, s.r.o.

Headquaters, warranty & service center: Električná 6471,
Trenčín 911 01, SK

CRN: 46259317
VAT ID: SK2023299729

Phone: +421(0) 917 700 098
E-mail: objednavky@insportline.sk
reklamacie@insportline.sk
servis@insportline.sk

Web: www.inSPORTline.sk

About shipping

