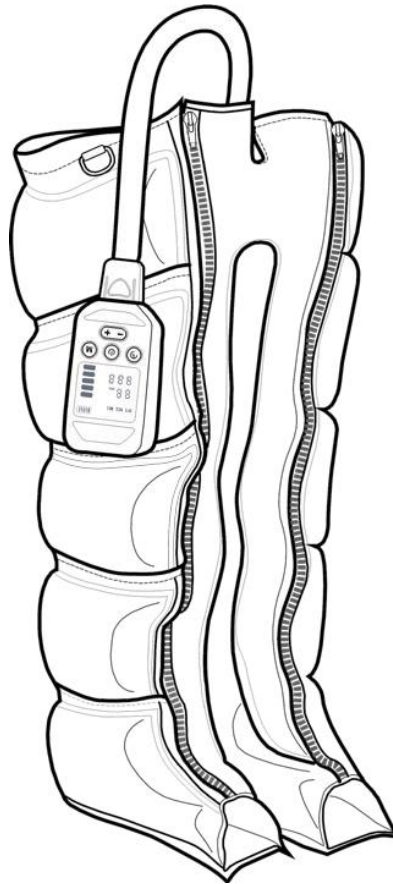




USER MANUAL – EN

IN 30014 Air Compression Leg Massager inSPORTline Complex



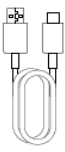

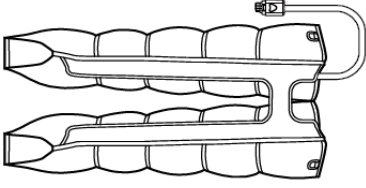


CONTENTS

- SAFETY INSTRUCTIONS..... 3
- PACKAGE CONTENTS 4
- PRODUCT DIMENSIONS 4
- USE 4
- COMPRESSOR UNIT CONNECTION 5
- CONTROLS..... 6
- TYPES OF MESSAGES..... 7
- CLEANING 7
- STORAGE 7
- FAQ 8
- ENVIRONMENT PROTECTION 8
- TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS..... 8


SAFETY INSTRUCTIONS

- Read the manual before use and keep it for future reference.
- Long-term continuous use can lead to muscle pain. Adjust the duration and intensity of use according to your body's reaction.
- Protect against impacts, drops, tears and similar physical damage. Keep away from fire, heat sources and flames.
- Protect against water and moisture. Store in a dry, shady and clean place. Do not store near heat sources or chemicals.
- Do not stand or perform other physically demanding activities during use.
- Charge only with a suitable cable and adapter.
- Do not charge unattended. Charge away from flammable or explosive substances.
- Do not cover the massager during use.
- The product is not suitable for pregnant women or people taking medications that could affect pressure, sensitivity, or judgment.
- The product is not suitable for people suffering from: hypertension, osteopenia, osteoporosis, diabetes, varicose veins, poor skin sensitivity, pressure sensitivity, recent injury or surgery, scoliosis or spinal deformity, people with heart disease or high blood pressure.
- Consult your doctor about use, especially if you suffer from the following medical complications or fall into the listed category: skin rash, open wounds, blisters, local tissue inflammation, bruising; deep vein thrombosis, osteomyelitis; fracture or myositis ossificans; high blood pressure; acute or severe heart, liver or kidney disease; neurological conditions leading to loss or change in sensitivity; bleeding disorders; recent surgery or injury at the application site; connective tissue disorders; chronic venous insufficiency; severe scoliosis or spinal deformity; you are taking a pacemaker, ICD or have had a heart attack; you are taking medications that thin the blood or alter sensitivity; you suffer from an allergy to the material of the product (nylon); you are pregnant or breastfeeding; you suffer from numbness; cancer or tumors; Epilepsy; cardiopathy
- Do not use the device on the face, neck, genitals.
- If you experience pain or discomfort, stop using the massager immediately.
- Do not use while charging.
- Charge regularly once every two months when the device is not in use.
- Do not wash in an automatic washing machine.
- The device is not a toy. Keep out of reach of children and pets.
- Charge the compression unit under supervision and away from flammable or explosive substances.
- Only charge the compression unit using the supplied charger.
- Do not dispose of the battery in a fire or leave it near heat sources.
- Do not dispose of the battery in water and protect against moisture.
- Do not repair or modify the device yourself.
- If there is an odor or liquid coming from the battery, stop operation immediately.
- Do not use the device in areas with high humidity and temperature.
- Suitable for home use only, not a substitute for health or medical care.

PACKAGE CONTENTS


				
Charging cable	Compression unit	Compression legs	Carrying bag	Manual

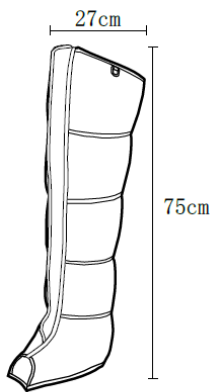PRODUCT DIMENSIONS




M

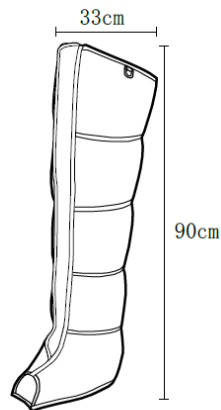Suitable for height 155-170 cm



L

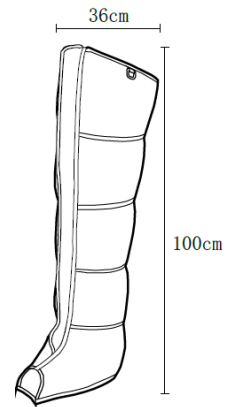Suitable for height 170-185 cm



XL

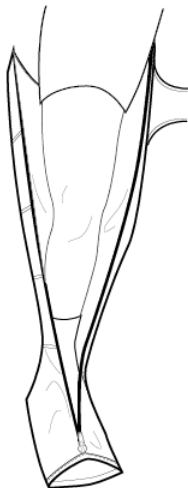
Suitable for height > 180 cm

USE

We recommend wearing shorts or thin pants during use. Sit or lie down comfortably. Your legs must be in line with your waist. Unzip the compression legs and place the legs separately in the leg. The leg should be in the center of the compression leg.


During use, the leg massager should not be wrinkled.

Do not use the device while standing or walking.




COMPRESSOR UNIT CONNECTION

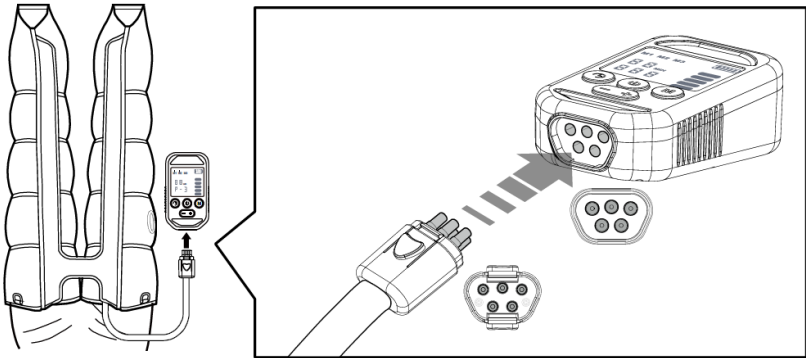Connect the inflation tubes to the air outlet of the compression unit. Make sure the connector is properly locked. You will hear a click.




You can attach the compression unit to the massager using a magnet.

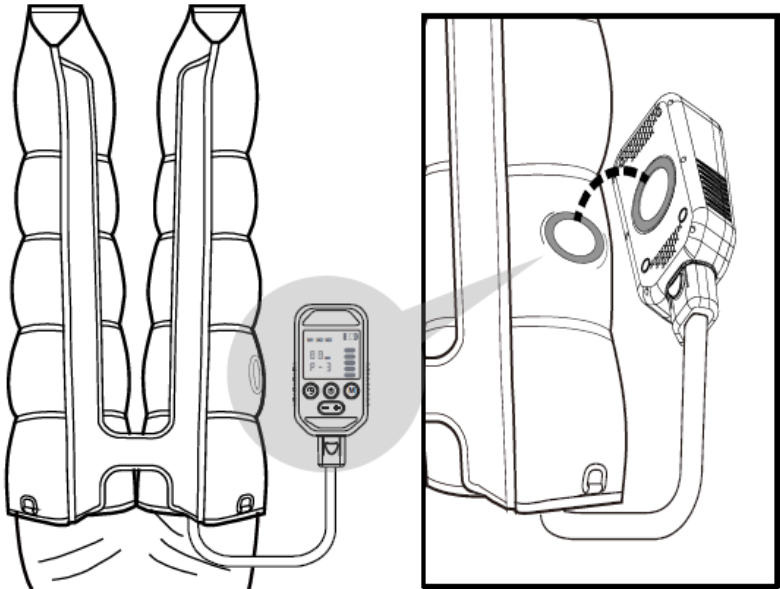Warning: do not place hands, fingers or other objects between the magnets, there is a risk of pinching.




Make sure the massager is closed and secured before starting. Do not use when the zipper is open.

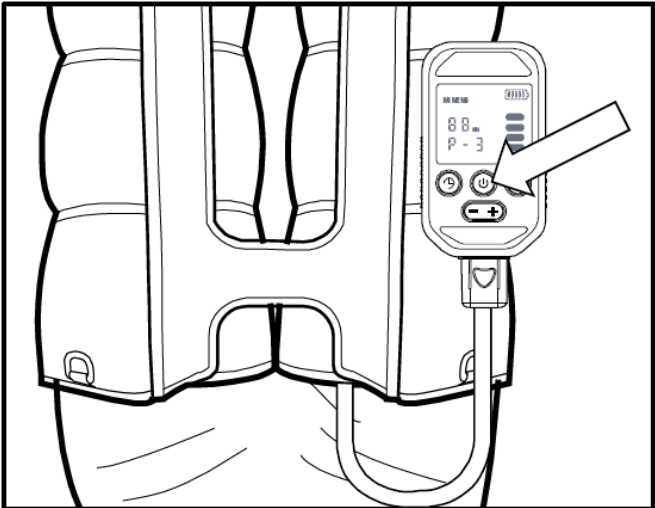
To start, press the power button for 1 second.



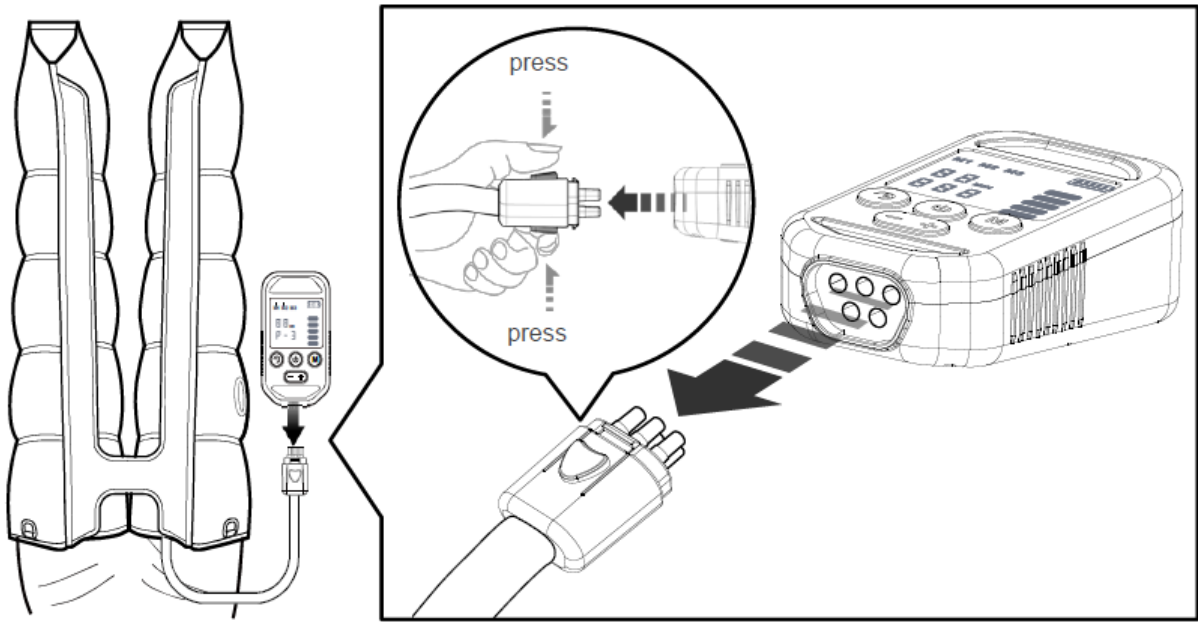To remove the compression unit, press the buttons on either side of the connector.



CONTROLS

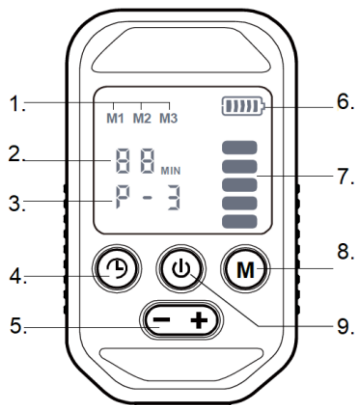


- | | |
|----|----------------------------------|
| 1. | Current Mode |
| 2. | Timer |
| 3. | Intensity |
| 4. | Timer button |
| 5. | Intensity button |
| 6. | Battery |
| 7. | Current leg compression position |
| 8. | Mode button |
| 9. | Power button |

The battery indicator shows 5 segments of the battery charge status. If the indicator flashes, the battery is low. The compression unit cannot be used while charging.

Button functions

Switches		Hold to start/stop the compression unit
Timer		Press to set the timer for 10, 20, 30 min.
Mode		Press to change the mode to M1, M2, M3. Hold to switch to manual mode, then press the button to adjust the compression position.
Intensity		Press to adjust the intensity P1 – P6. Press to decrease the intensity, press to increase it.

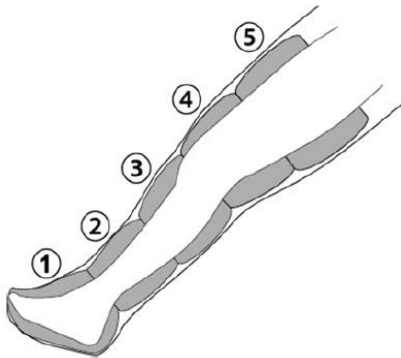TYPES OF MESSAGES



1. Sequential massage: the massage takes place gradually from the 1st to 5th chamber

2. Cyclic massage: cyclic massage of set chambers

3. Massage of a given point: setting the massage with a set chamber. Hold down the M button and then press to select a chamber.

MANUAL MODE

Hold down the M button to hear a beep.

In manual mode, you can select a massage with a specific chamber 1 – 5 or activate all chambers by pressing M.

Press to **■** decrease the intensity, press **+** to increase the intensity. P1 – P6

Using the timer button, you can then set the duration of use in the cycle 10, 20, 30 min.

In the event of a malfunction or emergency shutdown, you can detach the trouser legs from the compression unit.

CLEANING

- Spread the trousers on a flat surface. Use alcohol-based disinfectants for cleaning (do not use an oil-based agent). Do not use aggressive cleaning agents and solvents.
- Clean the outer parts with a soft cloth.
- You must dry the trousers properly before using them again.

Note: Do not immerse the massager in water or clean them under running water. Protect against the influence of aggressive cleaning agents and solvents.

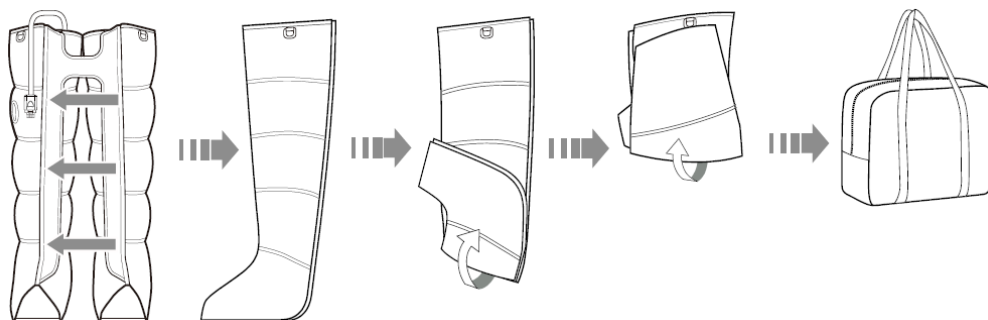
STORAGE

Spread on a flat, clean surface and let the trousers deflate.

Fold your right leg over your left. Align the edges.

From part of the feet, start folding the massager to one third and then again. Do not fold the massager too tightly, there is a risk of damaging the hoses and connections. Store the folded massager in a portable bag.

Make sure the trousers are clean and dry before storing. Store in a dry and shady place. Do not leave at places with high temperatures exceeding 45°C.



FAQ

Cannot be started.

Check the battery level. Make sure you press the power button correctly.

The display does not show the charge status while charging.

Make sure the connector is plugged in correctly.

Sudden shutdown during use.

Check the battery charge level. Shutdown to protect against overheating.

ERROR CODES



E-1 – Critically low battery

Charge the battery as soon as possible.



E-2 – air leak

Check the condition of the air hoses, their connections and all connections.

ENVIRONMENT PROTECTION

After the product lifespan expires or if the possible repairing is uneconomical, dispose of it according to the local laws and is environmentally friendly in the nearest scrapyard.

By proper disposal you will protect the environment and natural resources. Moreover, you can help protect human health. If you are not sure of the correct disposing, ask local authorities to avoid law violations or sanctions.

Don't put the batteries among house waste but hand them in to the recycling place.

TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS

General Conditions of Warranty and Definition of Terms

All Warranty Conditions stated here under determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Kutnohorská 531, Štěrboholy Retail Park, 109 00, Prague 10, ČR, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations, it depends whether the Buyer is the End Customer or not.

"The Buyer who is the End Customer" or simply the "End Customer" is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

"The Buyer who is not the End Customer" is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

Warranty Conditions

Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

Batteries

6-month battery warranty – we guarantee that battery's nominal capacity does not fall below 70% of its total capacity within 6 months of the product's sale.

The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damage
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints, wear of brake pads/blocks, chain, tires, cassette/multi wheel etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.

CZSEVEN SPORT s.r.o.

Registered Office: Kutnohorská 531, Štěřboholy Retail Park,
109 00, Prague 10, Czech Republic

Headquarters: Dělnická 957, Vítkov, 749 01Warranty &
Service: Čermenská 486, Vítkov 749 01

CRN: 26847264VAT ID:CZ26847264

Phone: +420 556 300 970

Email: eshop@insportline.cz
reklamace@insportline.cz
servis@insportline.cz

Web: www.inSPORTline.cz

SK**Stores inSPORTline SK, s.r.o.**

Headquaters, warranty & service center: Električná
6471, Trenčín 911 01, SK

CRN: 46259317VAT ID:SK2023299729

Phone: +421(0) 917 700 098E-mail:
objednavky@insportline.sk
reklamacie@insportline.sk
servis@insportline.sk

Web: www.inSPORTline.sk

About shipping