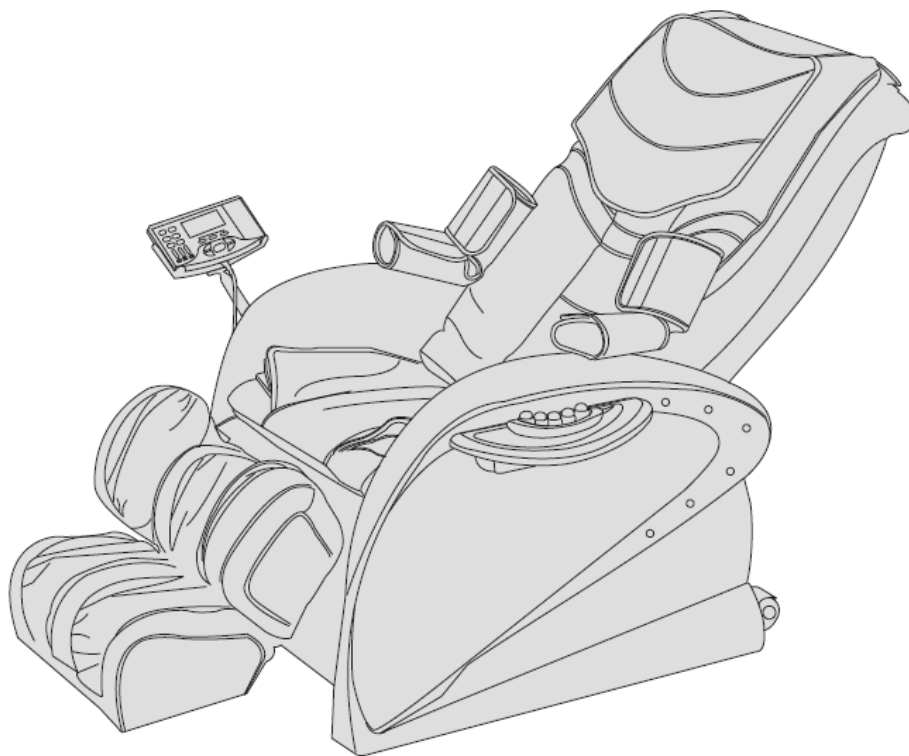




USER MANUAL – EN

IN 4477 Massage Chair inSPORTline SALLIERI



CONTENTS

SAFETY INSTRUCTIONS.....	3
PRODUCT SPECIFICATIONS.....	5
PARTS LIST	5
FUNCTIONS.....	6
Effects	7
CONTROL PANEL	7
Plug in	8
Adjusting backrest and footrest inclination	8
Start massage.....	8
Stop the massage, turn off the power	12
DOWNLOADING MP3 MUSIC	12
ADJUSTABLE PARTS	13
Padding and cushion assembly	13
Control panel and tray assembly	14
Armrest and footrest assembly	14
Armrest and footrest disassembly.....	15
Control panel holder assembly	15
Backrest electric cylinder assembly	17
PLACEMENT AND MOVING	17
Floor protection	17
Moving the chair.....	18
TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS.....	18

Thank you for purchasing this product. In order to use this product correctly, read this manual carefully before first use. It is important that you follow the safety precautions. Retain this manual for possible future use.

NOTE: We reserve the right to alter the design or the colour of the product without prior notice.

SAFETY INSTRUCTIONS

GENERAL INSTRUCTIONS

- Keep children away from the moving parts of the chair.
- Only use a power cable designed for this product and make sure it is grounded properly.
- To avoid an injury or damage, always plug the power cable out of the socket after use.
- Only use the chair as described in this manual. It is forbidden to use this chair in ways not listed in this manual.
- Do not use any accessories not recommended by the manufacturer.
- Do not use the chair outdoors.
- It is recommended to always use the chair 15 minutes at the time.
- Do not use the chair if the covering is damaged in any way.
- Do not use the chair with the air vent covered or if ventilation is obstructed in any way.
- Do not place any objects onto the chair.
- Make sure you do not fall asleep during the massage.
- Do not use the chair if you are under the influence of alcohol or if you are not feeling well.
- Do not use the chair at least one hour after your latest meal.
- Do not use a massage program that is too intense for you, you could sustain an injury.
- If the power cable is damaged, for safety reasons it should be replaced by the manufacturer, service technician or another qualified person.
- Maximum load: 150 kg

PLACEMENT

- Do not use the chair in high temperatures or in damp environment, for example in a bathroom.
- Do not use the chair directly after a significant change in temperature.
- Do not use the chair in dusty or otherwise unfavourable conditions.
- Keep a sufficient space around the chair and ensure proper ventilation.

RESTRICTIONS OF USE

- People with osteoporosis should not use this chair.
- People with heart disease or people with electronic implants, such as a pacemaker, should not use this chair.
- Do not use this chair if you have a fever.
- The chair should not be used by pregnant women or women on their period.
- Do not use this chair if you are injured or suffer from a skin condition.
- Children under the age of 12 should use the chair only under supervision.
- Physically or mentally handicapped people should use this chair only under the supervision of responsible and qualified person to ensure their safety.
- The chair should not be used by people in convalescence or people who are not feeling well.

- Do not use the chair if any part of your body is wet.
- Because of the heating surface of the chair, people with high sensitivity to heat should be careful while using this product.

SAFETY

- Make sure the outlet has the correct voltage.
- Don not insert the plug with wet hands.
- Water must not leak inside of the chair. The user may sustain an electric shock or the product could be damaged.
- Always make sure the cables are intact. Do not modify the electric circuit in any way.
- Do not use a wet cloth to clean the electronic parts, such as the power switch or the plug.
- If the power is cut off suddenly, do not approach the chair in case it's turned back on again. You could sustain an injury otherwise.
- If the chair stops to function properly, stop using it and contact the customer service centre.
- If you start to feel unwell, stop using the chair and consult your physician.

MAINTENANCE

- Maintenance should be carried out only by qualified service technicians. The user should not try to disassemble the chair or carry out maintenance by themselves.
- Don't forget to plug the device out of the socket after you finish using it.
- Do not use the chair if the plug or socket are damaged.
- Should the chair not be used for a longer period of time, curl up the cables and store the product in a dry and dust-free environment.
- Do not place the chair in an environment with high temperatures, in the vicinity of a fire or on direct sunlight.
- Clean this product with a dry cloth. Do not use alcohol-based cleaning products.
- The mechanical components of this product are specially made so that they do not need to be further maintained.
- Make sure the chair doesn't come into contact with any sharp objects.
- Do not roll or pull the chair across an uneven ground, it should always be lifted and carried.
- The chair is meant to be used intermittently. Do not use it for too long periods of time.


SOLUTIONS TO ORDINARY MALFUNCTIONS


- It is natural that the motor makes certain level of noise.
- If the chair doesn't work, make sure that the power cord is plugged in properly and that the power switch is on.
- If the maximum time of the massage is over, the system will automatically shut down. If the chair is turned on for too long, the temperature sensor safety will automatically shut it down. After that, the chair will need to be turned off for at least half an hour.

PRODUCT SPECIFICATIONS

Model	A27
Description	Massage chair
Input voltage	110 – 120 V – 60 Hz 220 – 240 V – 50 Hz 300 W
Safety class	Class I
Maximum time of usage	15 min
Used materials	PVC, PA, PU, leather, cloth, wood, steel and electronic parts
Weight limit	150 kg
Weight	112 kg
Dimensions after assembly	124 x 118 x 130 cm

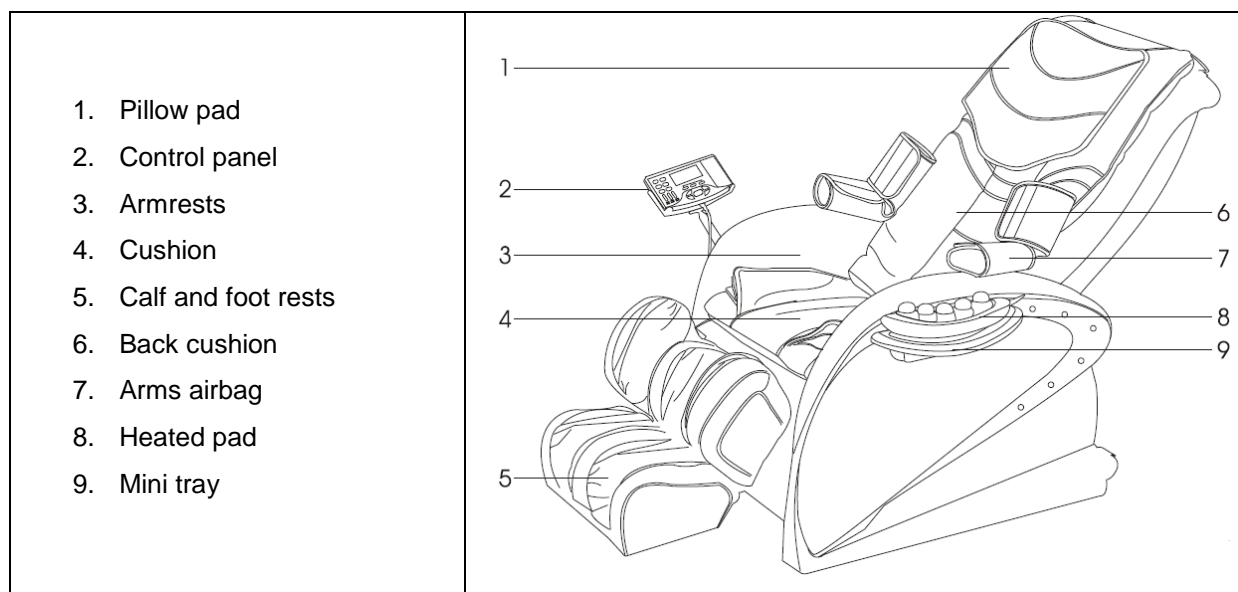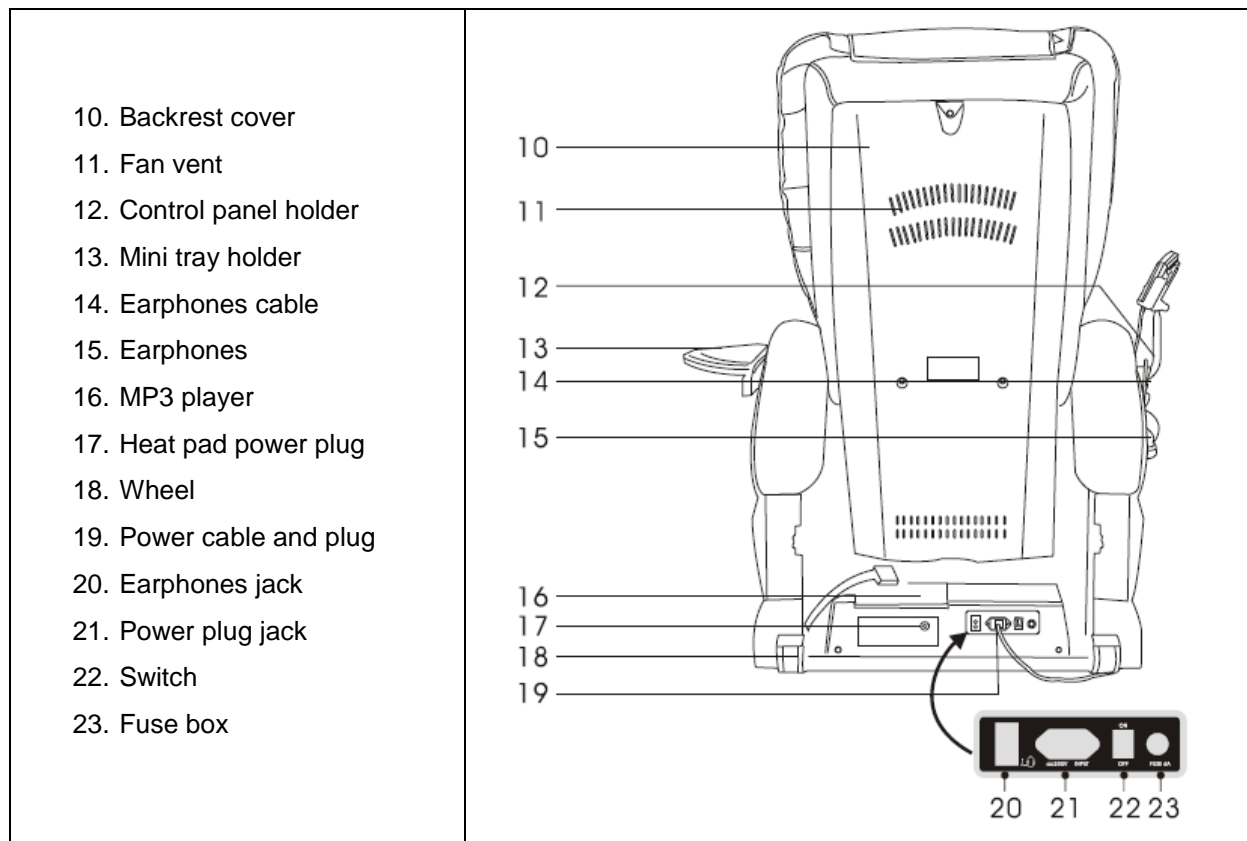PARTS LIST





FUNCTIONS

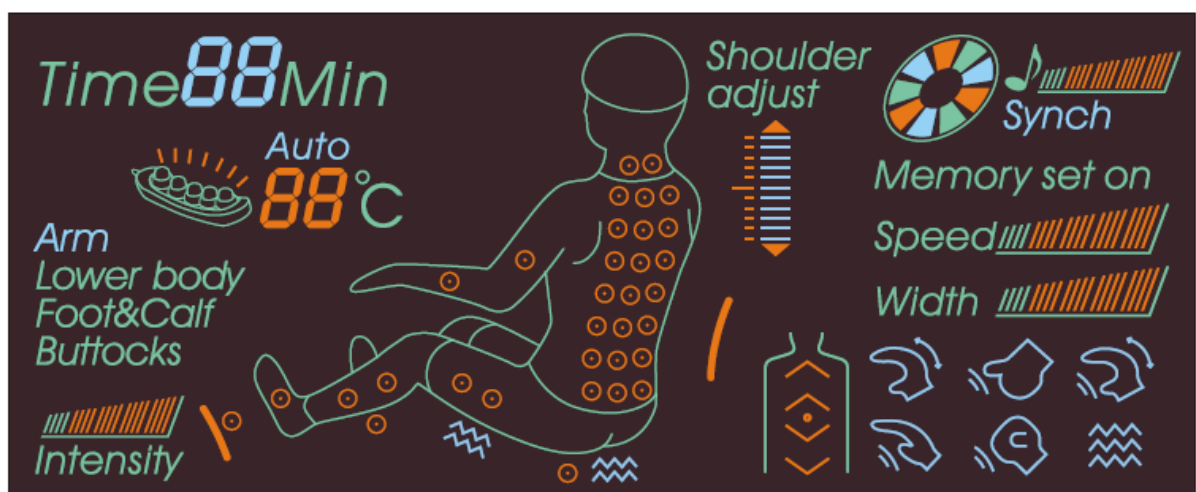
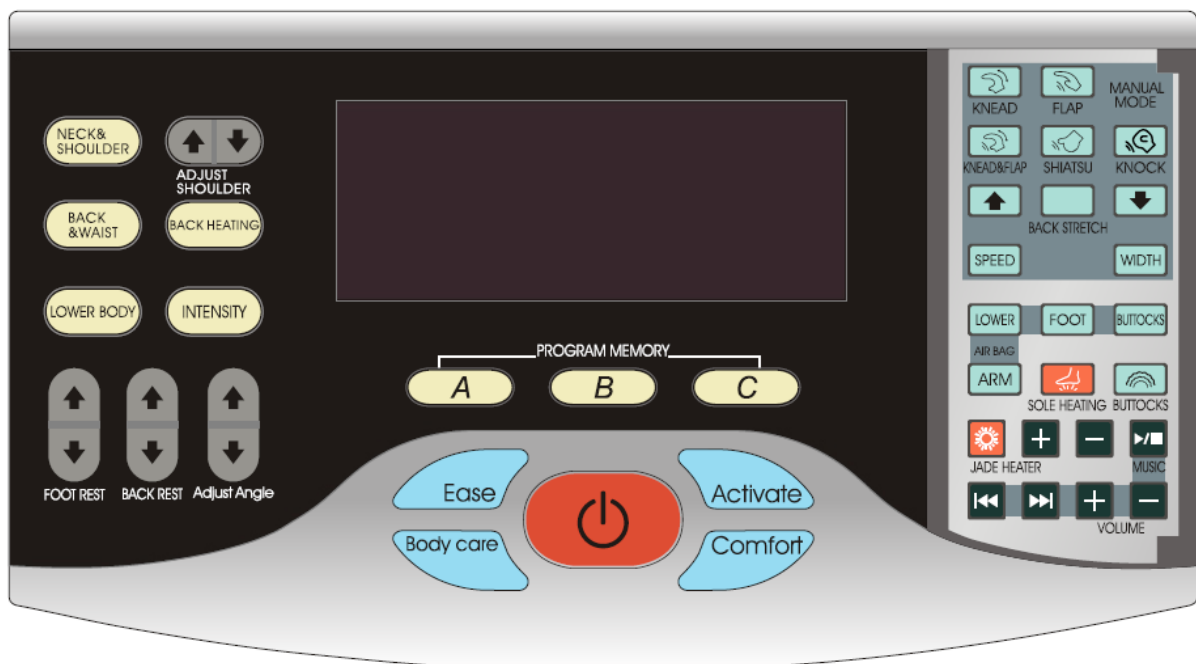
- The chair has a four-wheel quiet massage mechanism that moves the massage heads up and down in the backrest.
- **There are 5 massage emulation programs:** Shiatsu massage, kneading, flapping, knocking and synchronous kneading and flapping.
- **Pre-set automatic massage functions:** relax, pressure relieving, fitness and comfort.
- You can save three massages into the system's memory.
- Automatic massage of the upper body: shoulder and waist mode.
- There are three options for manual upper body massage – **overall, fixed area and partial area massage** – with **five levels of speed**. Under the flapping, knocking and shiatsu massage modes, you can set the distance between the massage heads in **three levels** – wide, middle and narrow.
- There are five airbags in the seat cushion for air pressure massage.
- The arms air massage is carried out through 4 inside airbags.
- **Vibration massage:** 3 levels of intensity in the buttocks area – strong, weak and off.
- **MP3 function:** The music can be played both during a massage and in the standby mode. You can adjust the volume or switch between songs.
- **Back heating massage:** The temperature of the cover can reach 40°C – 50°C; infrared heating.
- **Foot heating massage:** Infrared heating, maximum temperature 40°C.
- **Magnetic function:** In the calf and the sole area there are 8 magnets at 300 gauss strength.
- **Heating pad function:** Pre-set the temperature in the range of 40°C – 70°C. The temperature is adjusted by 5°C.

- **Chair inclination adjustment:** You can set the inclination of the back and foot rests.
- After switching the chair off, it will automatically return to its default position.
- The control panel has a VFD display.
- The chair's accessories include a MP3 player, headphones and a mini tray.

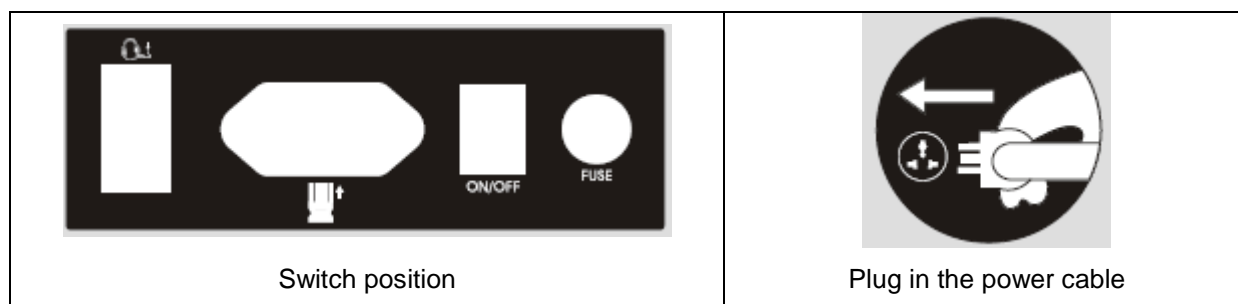
Effects

- Elimination of fatigue
- Blood circulation stimulation
- Relieving muscle fatigue
- Loosening tense muscles
- Relaxation
- Relieving muscle pain

CONTROL PANEL





Plug in





Adjusting backrest and footrest inclination

Do not adjust inclination during a massage.



Back rest adjustment:

Raise backrest		Press this button to raise the backrest slowly. Stop by releasing the button.
Lower backrest		Press this button to lower the backrest slowly. Stop by releasing the button.

Foot rest adjustment:

Raise footrest		Press this button to lower the footrest slowly. Stop by releasing the button.
Lower footrest		Press this button to lower the footrest slowly. Stop by releasing the button.



Whole frame adjustment:

Incline the whole frame backwards		Press this button to incline the frame slowly backwards. Stop by releasing the button.
Incline the whole frame forwards		Press this button to incline the frame slowly forwards. Stop by releasing the button.


Start massage

1. Press the red button to start the massage.
2. Detect your upper body massage point and then automatically detect your shoulder position. When you finish, you'll hear a chime. Now you can adjust your shoulder position. Press the massage button to start the massage. When adjusting the shoulder position, the A, B and C buttons will start flashing. With these buttons you'll select a pre-programmed massage mode. The buttons for automatic neck and shoulder massage, lower back massage and lower body massage will light up. If you press one of them, the current position will be saved.

3. You can adjust the massage heads positions. You have 5 positions to choose from.

	Move the massage head upwards.
	Move the massage head downwards.

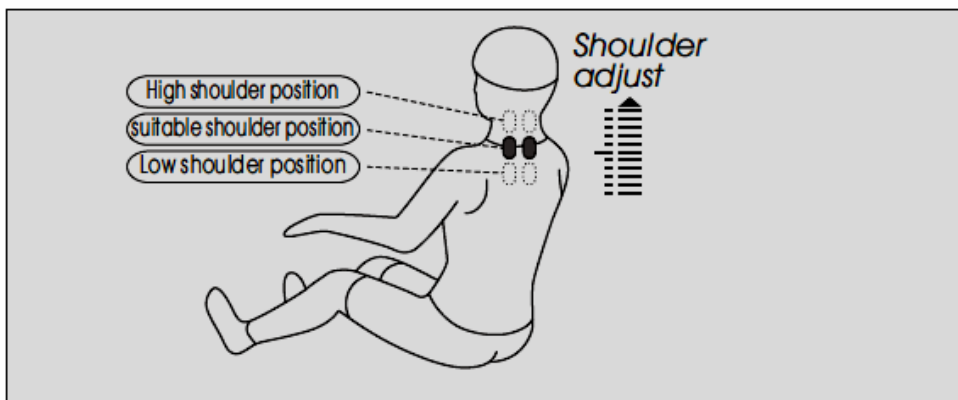The ideal massage position for your shoulders:










4. Press the message button to start the massage.

a. Automatic pre-set massage function:









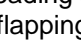







While setting the shoulder position, the A, B and C buttons will start flashing. If you press one of them, you'll begin a massage with saved settings. The memory function button will remain lit.



b. Automatic upper and lower body massage functions:


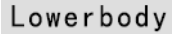



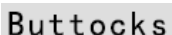
Name	Button	Message function	Effect	Display status
Automatic massage function		Comfort	Stimulates blood circulation	Auto 01
		Relax	Relaxes the whole body	Auto 02
		Fitness	Deep focused massage, relaxes aching muscles	Auto 03
		Pressure relieving	Relieves fatigue, makes the body active	Auto 04
Automatic massage		Neck and shoulder massage	Massage focused on the neck and shoulder area	Auto F1
		Back and waist massage	Massage focused	Auto F2
		Lower body massage, intensity can be adjusted	Massage focused on feet, calves, legs and buttocks	Lower body



c. Manual massage of the upper body:

Function	Button	Display status	Massage heads adjustment
Massage type	Knocking 		Speed, strength, width
	Shiatsu 		Speed, strength, width
	Kneading 		Speed, strength
	Flapping 		Speed, strength, width
	Kneading and flapping 		Speed, width
Speed		5 speed levels	
Width (range)		3 width levels	
Backrest stretch		Fixed position massage	
		Small area (partial) massage	
		Upper body (overall) massage	








Massage position can be adjusted with  a  buttons.







d. Manual lower body air pressure massage:











Function	Button	Description, display status	Options
Lower body massage		Air pressure massage of the lower body 	2 intensity levels
Feet massage		Air pressure massage of feet and calves 	2 intensity levels
Buttocks massage		Air pressure massage of buttocks 	2 intensity levels

Intensive massage		Intensity adjustment for lower body massage 	
----------------------	---	---	--

e. Other manual functions:


Function	Button	Description, status display	Options
Vibration		Vibration massage of buttocks 	2 intensity levels
Arm		Arm pressure massage	Massage of both arms
Backrest heating		Infrared heating at the back. The heat will become detectable by touch after 3 minutes. 	
Feet area heating		Infrared heating of the feet area. 	

Function	Button	Description, status display	Options
Heating pad		Heating pad switch 	
		Button to increase the temperature. The highest limit is 70 °C 	Increase by 5°C
		Button to decrease the temperature. The lowest limit is 40 °C 	Decrease by 5 °C

Function	Button	Description, display status	Options
Play music		Music on/off 	You can pause whenever you want
		Previous song 	
		Next song 	
		Volume up 	
		Volume down 	

Stop the massage, turn off the power

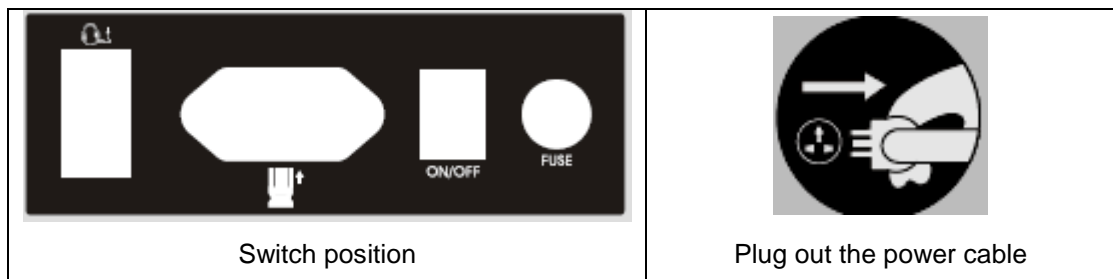1. Press the red power button to stop the massage. The backrest and footrest will automatically return to default position.
2. The memory buttons A, B and C will start flashing again. Press any of them to save the massage position before turning the chair off. The display will show "Memory Set On". Then press the red button to turn the chair off or wait 10 minutes for the chair to turn off automatically. The chair will revert to its default position after turning the power of.
3. Plug out the power cable:




DOWNLOADING MP3 MUSIC

1. Stop the massage and turn off the chair.
2. Pull the USB disk out of the MP3 player.
3. Insert the USB disk into a computer.
4. Download your favourite MP3-format music onto the USB disk. If you are unable to do so, format the disk.
5. After the download is complete, safely remove the device from the computer.
6. Insert the USB disk into MP3 player USB port.
7. Now you may turn the chair on again.


NOTE:

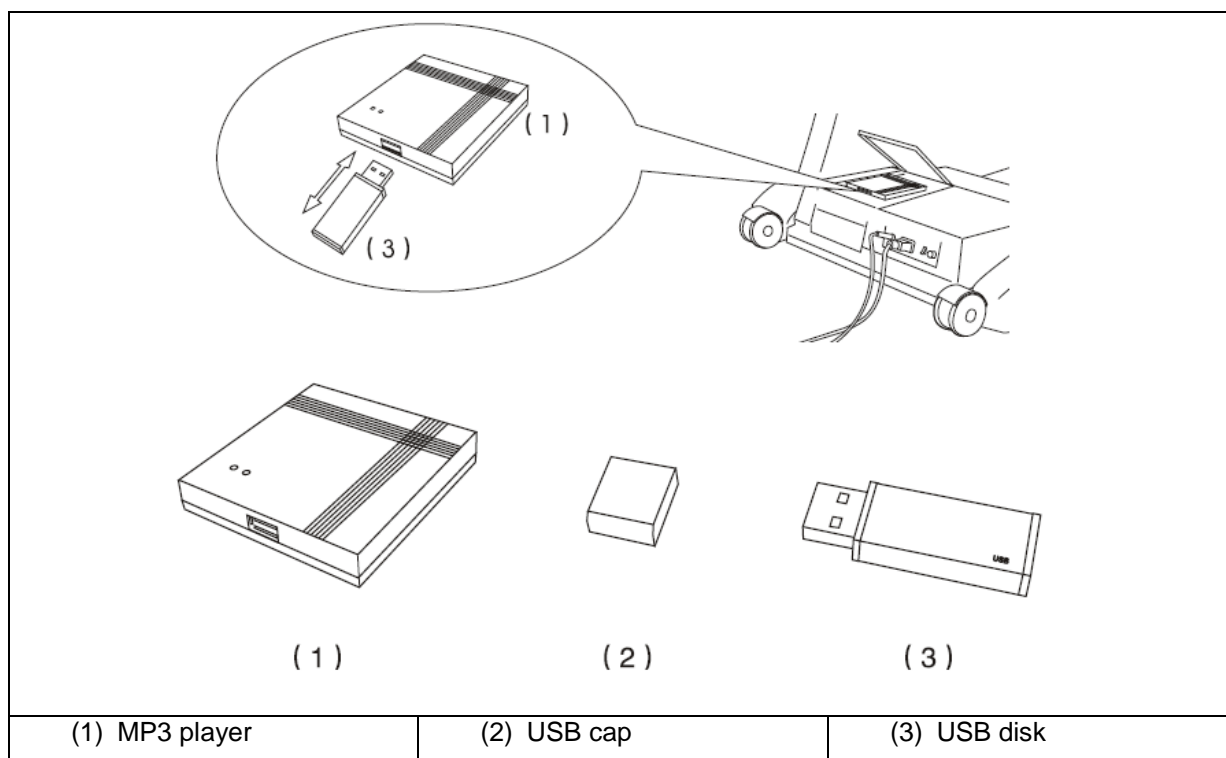1. If the MP3 player isn't able to play music from the USB disk, format the USB disk before downloading the music.
2. Only remove the USB disk, not the entire MP3 player.

**ADJUSTABLE PARTS****Padding and cushion assembly**


The padding can lessen the impact of the kneading in the neck and shoulder area. Attaching the padding is optional (it is recommended to use it).

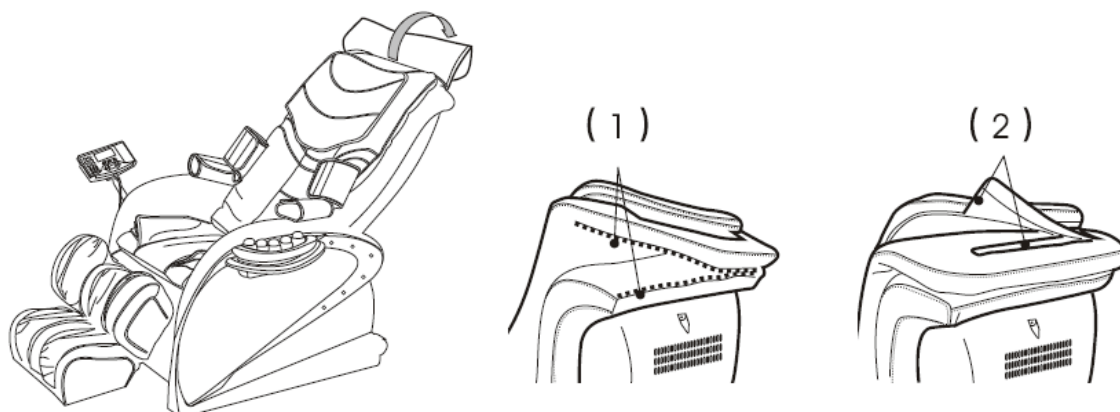Attach the padding according to the image – first, use a zipper to attach the upper part (1) and then attach the remaining part (2).




Control panel and tray assembly

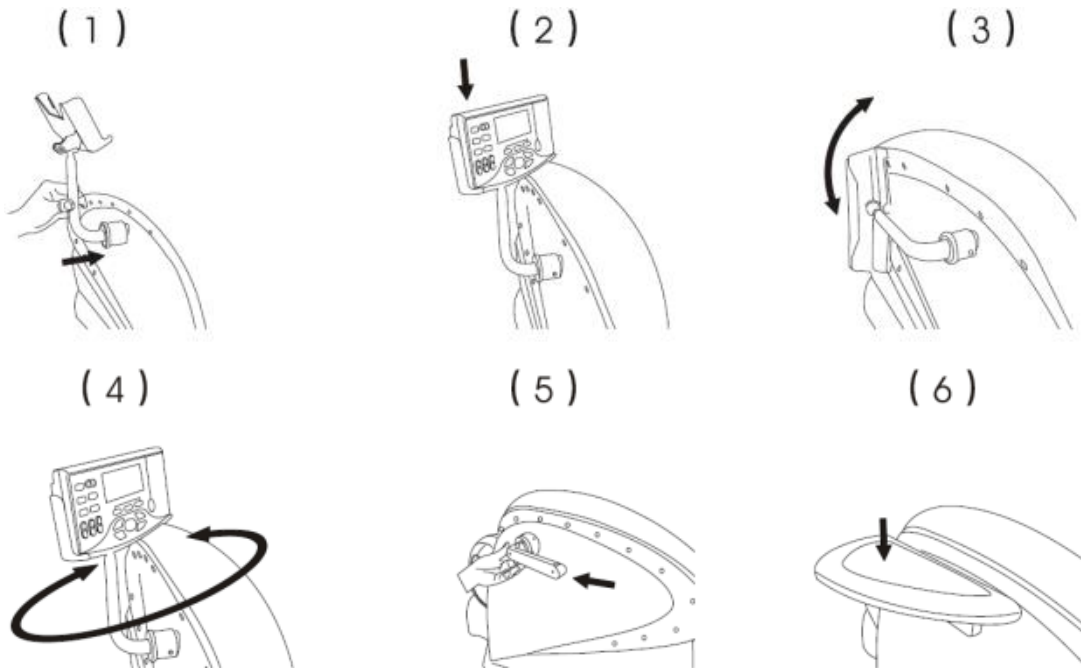Attach the control panel holder to the massage chair (1) and then insert the control panel itself (2). It should be easy to move the control panel both vertically (3) and horizontally (4). Then attach the tray holder (5) and insert the tray (6).




Armrest and footrest assembly

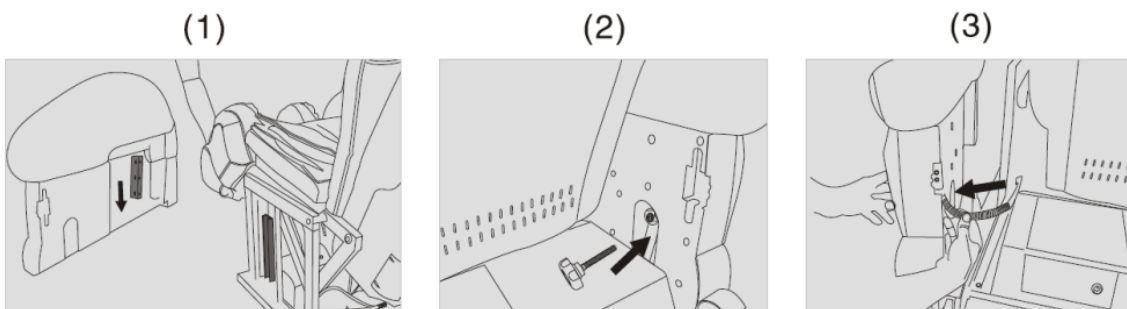1. One person should lift the calf rest, while another lifts the armrest. Let the armrest slide into the grooves on the frame (1).
2. Tighten the screw under the back of the armrest (2).
3. Affix the armrest air bag to the armrest and then cover the connection with a small cover board (3).







Armrest and footrest disassembly

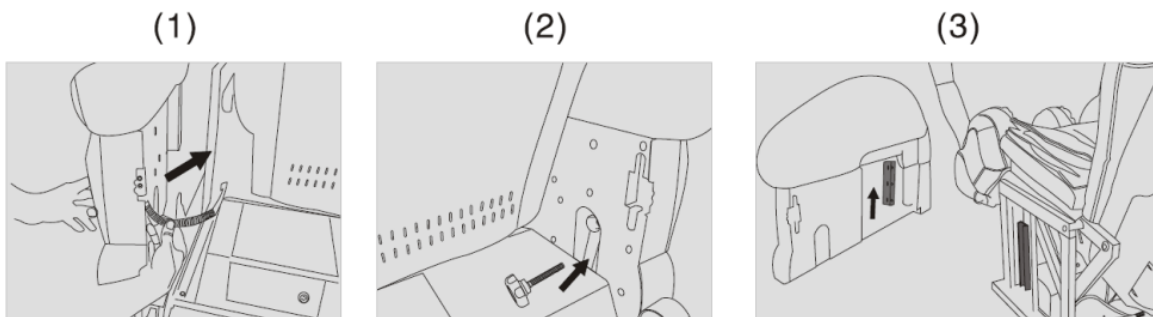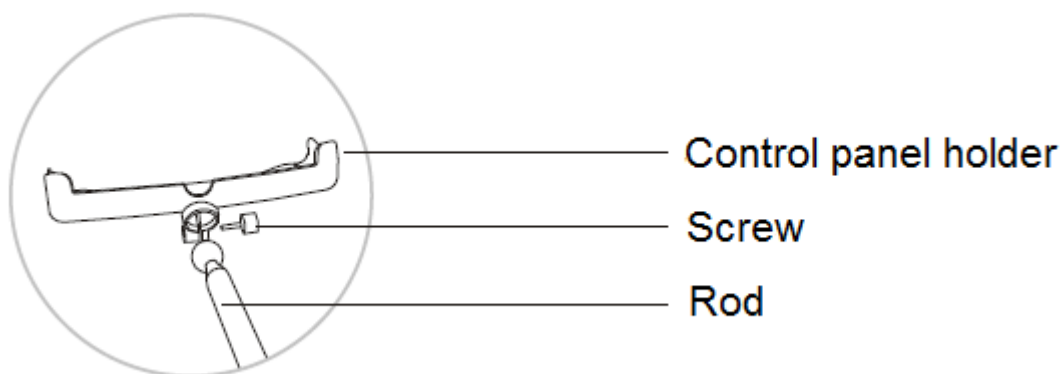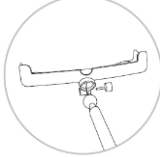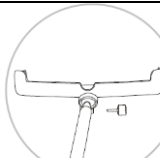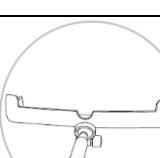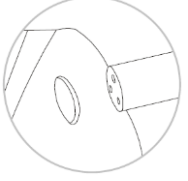
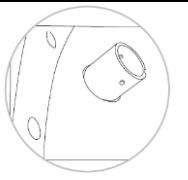

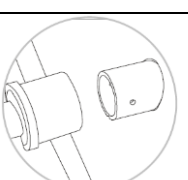
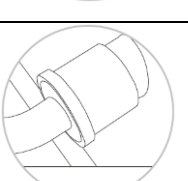
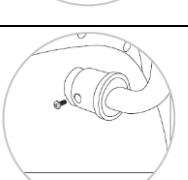
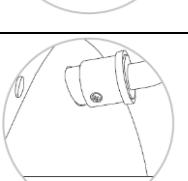
1. Remove the board cover and then detach the airbag connection (1).
2. Loosen the screw (2).
3. One person should lift up the footrest while another person lifts the armrest. Slide the armrests out of the frame (3).



Control panel holder assembly




	1. First, lay out all the parts and tools you'll need during the assembly.
	2. Insert the rod into the control panel holder.
	3. Tighten the screw clockwise.

	<p>4. Prepare the bushing.</p>
	<p>5. Insert the bushing into the hole on the right armrest so that the three screw holes on the bushing and armrest are aligned.</p>
	<p>6. Tighten the screws.</p>
	<p>7. Prepare the holder you assembled in step 3.</p>
	<p>8. Insert the holder into the bushing so that the screw holes on the holder and the bushing align.</p>
	<p>9. Prepare the screws.</p>
	<p>10. Tighten the screws.</p>

Backrest electric cylinder assembly

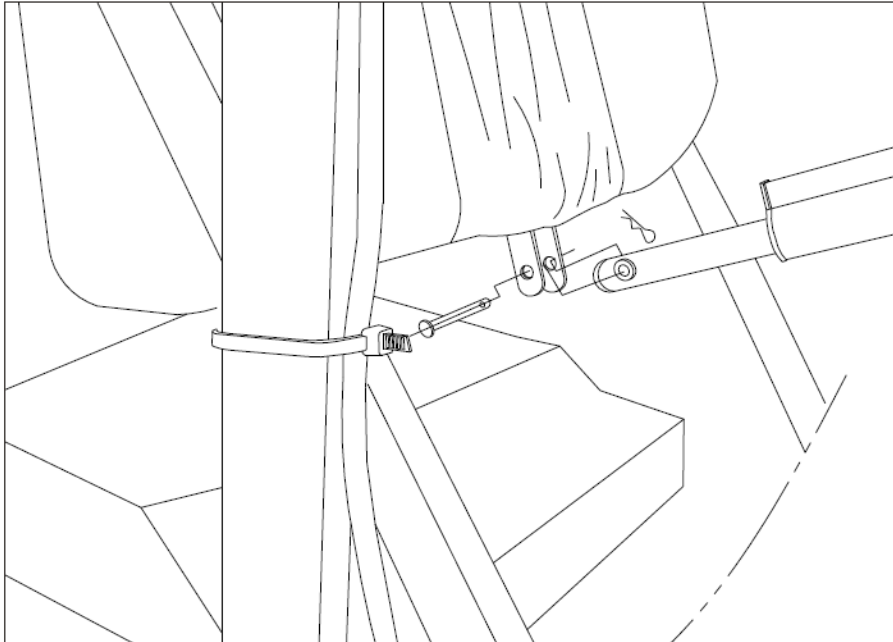1. Locate the support of the lower end of the backrest and place the electric cylinder into correct position.
2. Connect the support and the electric cylinder with a pin.
3. Secure the pin with a cotter.




PLACEMENT AND MOVING

Floor protection


If you place the massage chair on a wooden floor, the floor may get damaged after a period of time. Place the chair on a carpet or a protection pad.

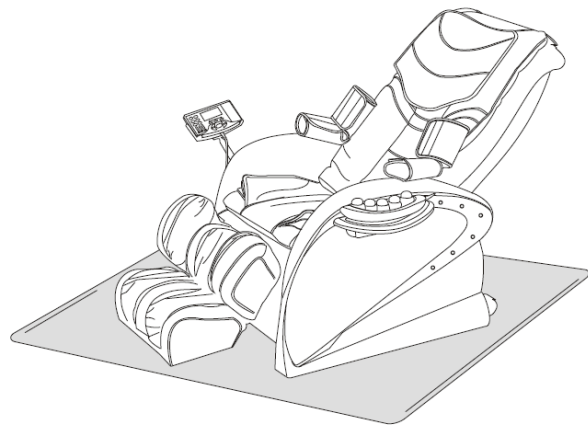WARNING: Do not move the chair on a wooden floor, uneven surface or in narrow spaces. In these cases, the chair should be lifted, by two people at least.



Moving the chair

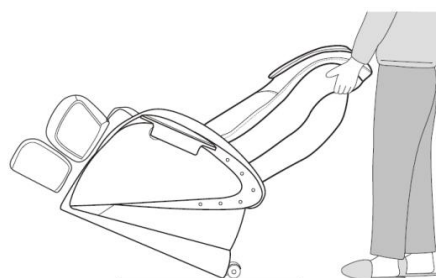Make sure the wires are not in the way before you tip the chair backwards so that it will stand on its wheels only. Then you can proceed to move the chair.



TERMS AND CONDITIONS OF WARRANTY, WARRANTY CLAIMS

General Conditions of Warranty and Definition of Terms

All Warranty Conditions stated hereunder determine Warranty Coverage and Warranty Claim Procedure. Conditions of Warranty and Warranty Claims are governed by Act No. 89/2012 Coll. Civil Code, and Act No. 634/1992 Coll., Consumer Protection, as amended, also in cases that are not specified by these Warranty rules.

The seller is SEVEN SPORT s.r.o. with its registered office in Borivojova Street 35/878, Prague 13000, Company Registration Number: 26847264, registered in the Trade Register at Regional Court in Prague, Section C, Insert No. 116888.

According to valid legal regulations it depends whether the Buyer is the End Customer or not.

“The Buyer who is the End Customer” or simply the “End Customer” is the legal entity that does not conclude and execute the Contract in order to run or promote his own trade or business activities.

“The Buyer who is not the End Customer” is a Businessman that buys Goods or uses services for the purpose of using the Goods or services for his own business activities. The Buyer conforms to the General Purchase Agreement and business conditions.

These Conditions of Warranty and Warranty Claims are an integral part of every Purchase Agreement made between the Seller and the Buyer. All Warranty Conditions are valid and binding, unless otherwise specified in the Purchase Agreement, in the Amendment to this Contract or in another written agreement.

Warranty Conditions

Warranty Period

The Seller provides the Buyer a 24 months Warranty for Goods Quality, unless otherwise specified in the Certificate of Warranty, Invoice, Bill of Delivery or other documents related to the Goods. The legal warranty period provided to the Consumer is not affected.

By the Warranty for Goods Quality, the Seller guarantees that the delivered Goods shall be, for a certain period of time, suitable for regular or contracted use, and that the Goods shall maintain its regular or contracted features.

The Warranty does not cover defects resulting from (if applicable):

- User's fault, i.e. product damage caused by unqualified repair work, improper assembly, insufficient insertion of seat post into frame, insufficient tightening of pedals and cranks
- Improper maintenance
- Mechanical damages
- Regular use (e.g. wearing out of rubber and plastic parts, moving mechanisms, joints etc.)
- Unavoidable event, natural disaster
- Adjustments made by unqualified person
- Improper maintenance, improper placement, damages caused by low or high temperature, water, inappropriate pressure, shocks, intentional changes in design or construction etc.

Warranty Claim Procedure

The Buyer is obliged to check the Goods delivered by the Seller immediately after taking the responsibility for the Goods and its damages, i.e. immediately after its delivery. The Buyer must check the Goods so that he discovers all the defects that can be discovered by such check.

When making a Warranty Claim the Buyer is obliged, on request of the Seller, to prove the purchase and validity of the claim by the Invoice or Bill of Delivery that includes the product's serial number, or eventually by the documents without the serial number. If the Buyer does not prove the validity of the Warranty Claim by these documents, the Seller has the right to reject the Warranty Claim.

If the Buyer gives notice of a defect that is not covered by the Warranty (e.g. in the case that the Warranty Conditions were not fulfilled or in the case of reporting the defect by mistake etc.), the Seller is eligible to require a compensation for all the costs arising from the repair. The cost shall be calculated according to the valid price list of services and transport costs.

If the Seller finds out (by testing) that the product is not damaged, the Warranty Claim is not accepted. The Seller reserves the right to claim a compensation for costs arising from the false Warranty Claim.

In case the Buyer makes a claim about the Goods that is legally covered by the Warranty provided by the Seller, the Seller shall fix the reported defects by means of repair or by the exchange of the damaged part or product for a new one. Based on the agreement of the Buyer, the Seller has the right to exchange the defected Goods for a fully compatible Goods of the same or better technical characteristics. The Seller is entitled to choose the form of the Warranty Claim Procedures described in this paragraph.

The Seller shall settle the Warranty Claim within 30 days after the delivery of the defective Goods, unless a longer period has been agreed upon. The day when the repaired or exchanged Goods is handed over to the Buyer is considered to be the day of the Warranty Claim settlement. When the Seller is not able to settle the Warranty Claim within the agreed period due to the specific nature of the Goods defect, he and the Buyer shall make an agreement about an alternative solution. In case such agreement is not made, the Seller is obliged to provide the Buyer with a financial compensation in the form of a refund.



SEVEN SPORT s.r.o.

Registered Office:	Borivojova 35/878, 130 00 Praha 3, Czech Republic
Headquarters:	Delnicka 957, 749 01 Vitkov, Czech Republic
Warranty & Service Centre:	Cermenska 486, 749 01 Vitkov, Czech Republic
CRN:	26847264
VAT ID:	CZ26847264
Phone:	+420 556 300 970
E-mail:	eshop@insportline.cz reklamace@insportline.cz servis@insportline.cz
Web:	www.insportline.cz

SK

INSPORTLINE s.r.o.

Headquarters, Warranty & Service centre: Elektrická 6471, 911 01 Trenčín, Slovakia

CRN: 36311723

VAT ID: SK2020177082

Phone: +421(0)326 526 701

E-mail: objednavky@insportline.sk

reklamacie@insportline.sk

servis@insportline.sk

Web: www.insportline.sk

Date of Sale:

Stamp and Signature of Seller: